

CM200 Series

K-12 EDUCATIONAL COMPOUND MICROSCOPE

This versatile and value priced microscope is feature packed for use at all levels of K-12 education. Designed for many years of rugged use, the CM200 series is equipped with three DIN Standard objectives and Wide Field 10x eyepiece which ensure that every slide renders the best resolution possible. Available with either traditional tungsten illumination or the latest in cordless **Green Friendly** LED technology.

General features:

Each CM200 is die cast in aluminum and designed to withstands the constant rigors of the classroom. Separate coarse & fine focus controls for ease of use with a slip clutch to protect gears from eager students.

Optics:

High quality DIN standard optics in three powers with a retractable 40x for safety.

Illumination:

The CM200 is equipped in either the industry standard 20 Watt Tungsten bulb with dimmer switch or **Green Friendly** cordless LED with variable light intensity.

Quality Control:

Each CM200 is rigorously tested at the factory which is controlled by ISO9001 and ISO14001 system.

Service After Sales:

Complete inventory of parts and convenient repair service ensure that your investment is well taken care of.






Features At A Glance

- ◇ **W10X Eyepiece w/ wire pointer**
- ◇ **Locked on eyepiece, slide clips**
- ◇ **All glass lenses**
- ◇ **Cord wrap on corded model**
- ◇ **Slip clutch gear protector**
- ◇ **Plastic Iris diaphragm protector**
- ◇ **Built to last**
- ◇ **Green friendly LED illumination available with long life bulb**
- ◇ **Cool white light**

Nexcope

CM200 Series

Specifications

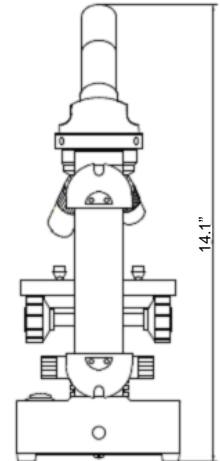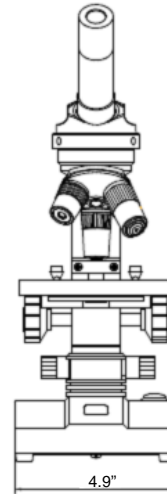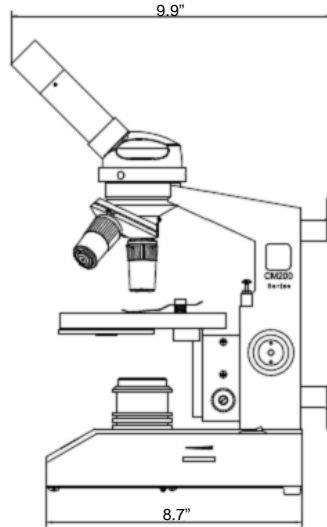

LENGTH	WIDTH	HEIGHT	WEIGHT
7.8"	4.9"	14.1"	6.75 lbs

MODEL	HEAD	EYEPIECE	OBJECTIVE	CONDENSER	ILLUMINATOR
CM201	Monocular	W10XD, 18mm	4XD, 10XD, 40XRD	Built-in NA 0.65 w/ disc diaphragm	120V, 20W w/ dimmer
CM202	Monocular	W10XD, 18mm	4XD, 10XD, 40XRD	Built-in NA 0.65 w/ iris diaphragm	Cordless LED w/ dimmer

Accessories for CM200 Series:

OBJECTIVES

CMA321 4XD Achromatic
 CMA322 10XD Achromatic
 CMA323 40XRD Achromatic


EYEPIECES

CMA275 W10XD, 18mm w/wire pointer

MISCEL:


CMA335 UL/CE approved power adapter/charger

ISO9001:2000

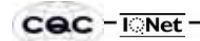

00109Q11543R2M/3

ISO14001:2004

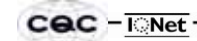

00107E20511R0M/330



Design Changes:
Nexcope
 reserves the right to make changes in instrument design in accordance with changes in technology without notice or obligation.

Authorized Dealer:

Ample Scientific L.L.C.
 4220 Steve Reynolds Boulevard, Suite 9, Norcross, GA 30093
 Phone: (678) 380-0369
 E-mail: sales@amplescientific.com
 Http://www.amplescientific.com